



# BOUGIE ASSISTED INTUBATION

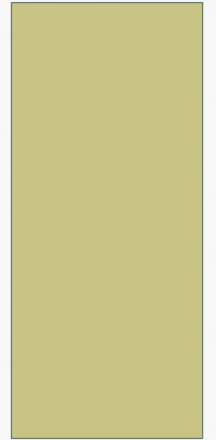


# BOUGIE ASSISTED

- We all know a bougie can make things easier when it comes to intubation
- Using it as a stylet with the ETT can be challenging to control
  - Bougie wraps around and is loose beyond the distal end of the ETT
  - Bougie can rotate making it harder to control

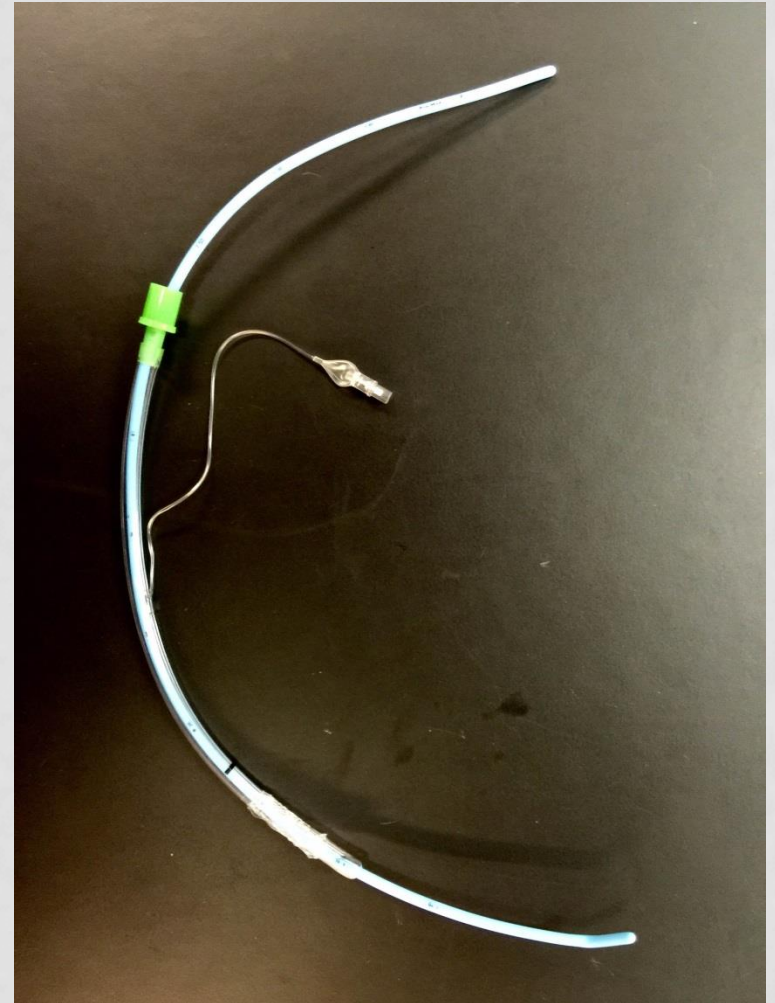


# BOUGIE ASSISTED

- How can we make this easier to control and use?
- Couple options
- 1:
  - Place the ETT over the bougie until the end of the ETT is at the coude tip of the bougie
  - Wrap the opposite or long end of the bougie back into the inlet of the ETT

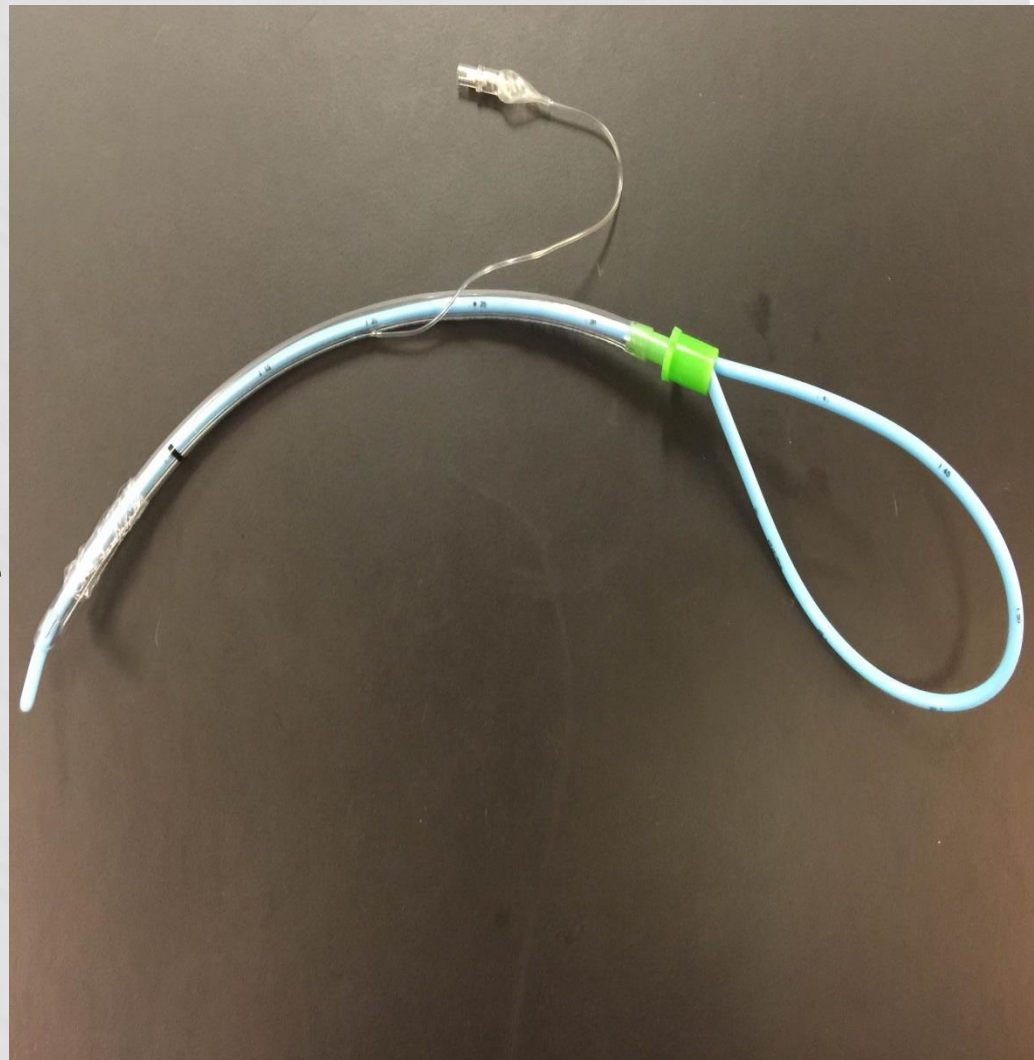

# BOUGIE ASSISTED

- 2:
  - Same setup as the last with ETT placed over the bougie down to the coude tip
  - Pop the 15mm connection apart from the ETT (as it commonly likes to anyways)
  - Place the long end of the bougie into the actual lumen of the ETT

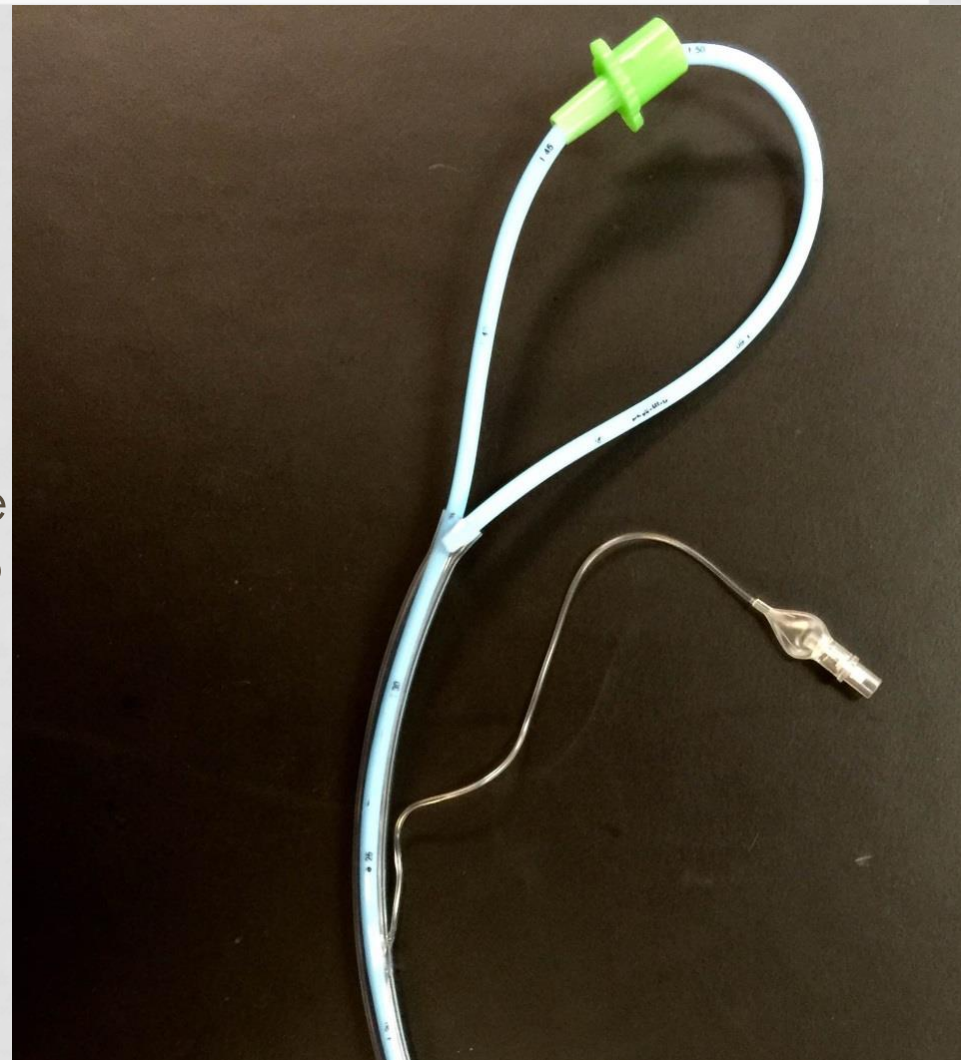

# BOUGIE ASSISTED

- These options give you better control of the bougie and prevent it from rolling over
- Also you can use it as you do without the ETT
  - Use a trigger method to raise the tip of the bougie